



Introducing the SureTemp[®] Plus 692 Electronic Thermometer

Speed and accuracy you require.
Safety and security you demand.

WelchAllyn[®]

Advancing Frontline Care[™]



The SureTemp® Plus 692 – with the speed and accuracy you require, and the safety and security you demand

- **Easy to use, intuitive operation**
- **Removable probe well reduces the risk of cross-contamination**
- **Improved accuracy for even more reliable readings**
- **Innovative security features that cut down on the chance of theft**

When we set out to improve the Welch Allyn SureTemp, we knew we had our work cut out for us. After all, when you start with the thermometer that helped set the standard for speed and accuracy, the expectations run pretty high.

So how did we do? See for yourself by taking a look at the innovative SureTemp Plus 692 – easier, safer and more accurate than the original, and still built to withstand constant use in demanding environments.

Easy: An easy-to-read LCD display with flashing icons that clearly communicate in any language, and an intuitive user interface and ergonomic shape make using the SureTemp Plus 692 as simple as possible.

Safe: Patient safety is improved thanks to an interchangeable probe well that reduces the risk of cross-contamination.

Fast and Accurate: In addition to improved accuracy, the SureTemp 692 is one of the fastest thermometers on the market, providing approximately 4- to 6-second oral, 10-second paediatric axillary, 15-second adult axillary, and 10-second rectal measurements.

Secure: To top it all off, our new security features make sure the unit stays where it is supposed to be. A location ID on the display and other innovative features help prevent theft.

ID Field for electronic labelling of thermometer



Battery life indicator

User selectable icon and/or words for identification of thermometer mode

Interchangeable probe with push-button probe cover ejection

Removable probe well for easy cleaning

Convenient storage for 25 probe covers


C/F button for toggling between Fahrenheit and Celsius

Large, easy-to-read, backlit LCD

Mode button for selecting oral, paediatric axillary, or adult axillary

Button for timing pulse and respiration

Ergonomic design with rubberised grips

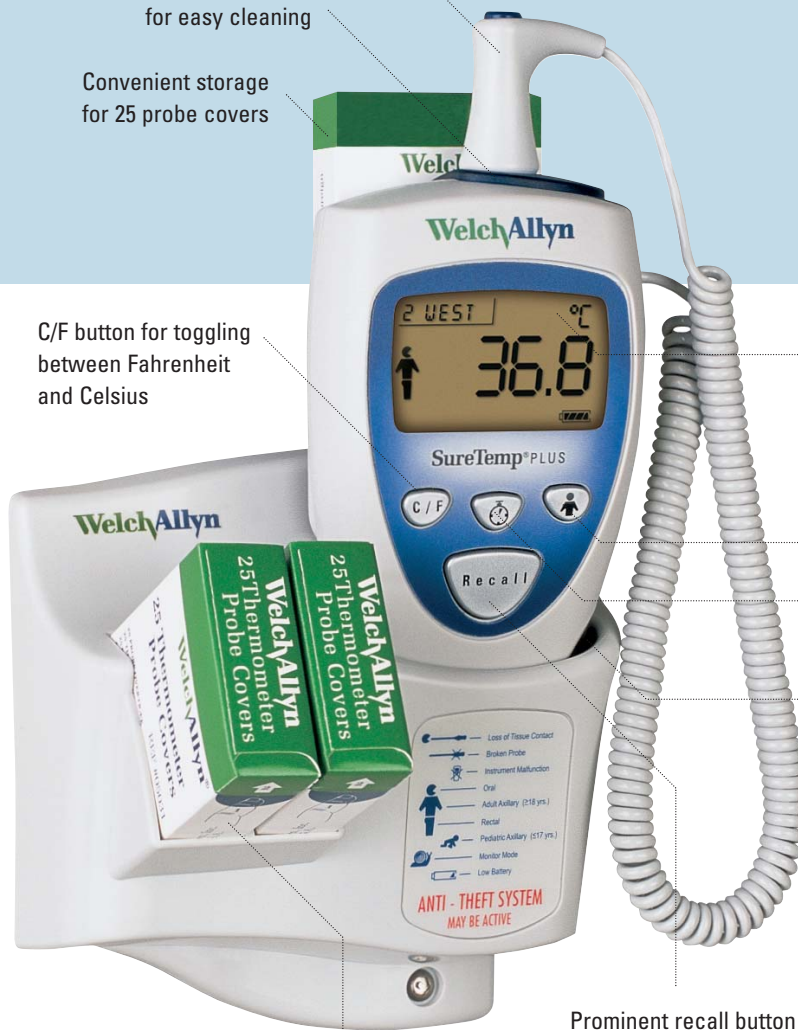

Prominent recall button for display of last temperature taken

Storage space for 50 additional probe covers in the wall mount (Wall mount standard with the 692 model)

Interchangeable rectal probe and well reduces risk of cross-contamination



Optional rolling stand for security and portability

The SureTemp® Plus 692

WelchAllyn®

Advancing Frontline Care™

The SureTemp® Plus 692 –
 setting the standard of
 care even higher



Welch Allyn U.K. Ltd.
 Cublington Road, Aston Abbots
 Buckinghamshire HP22 4ND, England
 Tel: 0 207 365 6780 Fax: 0 207 365 9694

Welch Allyn European Customer Service
 Navan Business Park, Dublin Road
 Navan, County Meath, Ireland
 Tel: +353 46 906 7790 Fax: +353 46 906 7754

Welch Allyn PTY. LTD, S. Africa
 394 Surrey Ave, Ferndale, Randburg 2194 South Africa
 P O Box 809, Fontainebleau, South Africa 2032
 Tel: 0 11 777 7555 Fax: 0 11 777 7556

SureTemp® Plus 692 specifications

Calibration Accuracy.....	±0.1 °C (Monitor Mode)
Temperature Predict Time:	
Oral.....	Approx. 4 – 6 secs
Adult Axillary (18 years and older).....	Approx. 12 – 15 secs
Paediatric Axillary (17 years and younger).....	Approx. 10 – 13 secs
Rectal.....	Approx. 10 – 13 secs
Unit Dimensions.....	215 mm x 81 mm x 62 mm
Weight.....	357 grams
Ambient Temperature Range.....	10 °C – 40.0 °C
Patient Temperature Range.....	26.7 °C – 43.3 °C
Power Source.....	(3)1.5Vdc AA batteries
Battery Life.....	Approx. 5000 temperatures
Warranty.....	3 year limited instrument, 1 year limited probe, and 90 day limited probe well

Ordering information

Thermometers

- 01692-400** SureTemp Plus Electronic Thermometer (Model 692), Wall Mount, Security System with ID Location Field, and Oral Probe with Oral Probe Well, 3 Year Limited Warranty
- 01692-401** SureTemp Plus Electronic Thermometer (Model 692), Wall Mount, Security System with ID Location Field, and Rectal Probe with Rectal Probe Well, 3 Year Limited Warranty
- 01692-500** SureTemp Plus Electronic Thermometer (Model 692), One Per Room, Wall Mount, 2.74 m Oral Probe with Oral Probe Well, 3 Year Limited Warranty (*must be locked to the wall*)
- 01692-501** SureTemp Plus Electronic Thermometer (Model 692), One Per Room, Wall Mount, 2.74 m Rectal Probe with Rectal Probe Well, 3 Year Limited Warranty (*must be locked to the wall*)

Accessories

- 21326-0000** Wall Holder
- 21330-0000** Holder for Spare Probe & Well
- 06138-000** Calibration Key Assy, 690/692
- 02892-000** Probe & Well Kit 1.22 m Rectal
- 02892-100** Probe & Well Kit 2.74 m Rectal
- 02893-000** Probe & Well Kit 1.22 m Oral
- 02893-100** Probe & Well Kit 2.74 m Oral
- 02891-0000** Probe Well, Blue
- 02891-1000** Probe Well, Red
- 08273-000** Thermometer Stand with Holder for Spare Probe & Well
- 21333-0000** Bracket, M600 Stand Upgrade, with Holder for Spare Probe & Well
- 98025-0000** Extended Warranty M692

Replacement probe covers

- 05031-101** Case of 1,000 probe covers

Call Welch Allyn to order the SureTemp® Plus 692
 4341 State Street Road, PO Box 220,
 Skaneateles Falls, NY 13153-0220, USA
 Tel: 800 535 6663 • Fax: 315 685 3361
www.welchallyn.com

WelchAllyn®

Advancing Frontline Care™